

ABIN1 (mouse; full length), pAb

Alternate Names: A20-binding inhibitor of NF-kappa-B activation 1, HIV-1 Nef-interacting protein, Nef-associated factor 1, Nip40-1, Virion-associated nuclear shuttling protein, TNFAIP3-interacting protein 1

Cat. No. 68-0001-100
Lot. No. 30238

Quantity: 100 µg
Storage: -20°C

FOR RESEARCH USE ONLY

NOT FOR USE IN HUMANS

CERTIFICATE OF ANALYSIS

Page 1 of 2

This antibody was developed and validated by the Medical Research Council Protein Phosphorylation and Ubiquitylation Unit (University of Dundee, Dundee, UK).

Background

ABIN1 contains a polyubiquitin-binding domain homologous to that present in nuclear factor κB (NF-κB) essential modulator (NEMO), a component of the inhibitor of NF-κB (IκB) kinase (IKK) complex (Nanda *et al.*, 2011). The interaction of ABIN1 with polyubiquitin is required to limit the activation of the TLR (Toll-like receptor) signalling pathway and prevent autoimmunity (Nanda *et al.*, 2011). ABIN1 has been implicated in autoimmune disorders such as systemic lupus erythematosus (SLE), diabetes mellitus and psoriasis (Gao *et al.*, 2011; Zhou *et al.*, 2011) and has the ability to inhibit TNF-induced apoptosis (Verstrepena *et al.*, 2009). The TLR pathway activates diverse signal-transduction pathways that control proinflammatory effector mechanisms, such as production of cytokines, growth factors and cytotoxic mediators (Zhou *et al.*, 2011).

Antibody Production:

Anti-ABIN1 (mouse) polyclonal antibody was raised in sheep against ABIN1 (residues 1-560 of mouse ABIN1). The antibodies were purified by the Medical Research Council Protein Phosphorylation and Ubiquitylation Unit (MRC-PPU, University of Dundee, Dundee, U.K.) by affinity purification of the anti-ABIN1 pAbs from the sheep serum using a GST-tagged antigen-agarose column. Anti-ABIN1 (mouse) pAb was sourced by Ubiquigent directly from the MRC-PPU.

Continued on page 2

Physical Characteristics

Quantity: 100 µg

Formulation: phosphate-buffered saline

Concentration: to be provided on shipping

Specificity: detects ABIN1 at ~73 kDa

Source: sheep polyclonal antibody

Reactivity: mouse; other species not tested

Immunogen: mouse ABIN1 (residues 1-560)


Stability/Storage: 12 months at -20°C; aliquot as required

Purification: affinity-purified using immobilized immunogen

Research Applications and Quality Assurance

Western Immunoblotting:
Use 0.3 µg/ml

Immunoprecipitation:
Use 0.3 µg/mg of cell extract

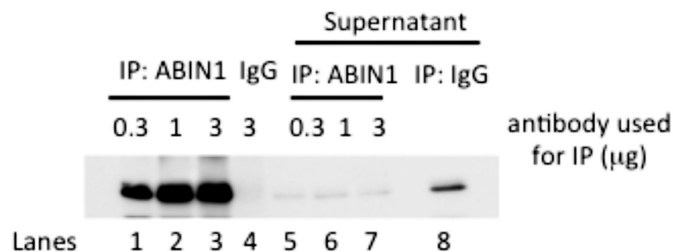

Immunoprecipitation Assay:

ABIN1 was immunoprecipitated from Raw 264.7 total cell extracts (1 mg) using various amounts of anti-ABIN1 antibody (Cat# 68-0001-100) or pre-immune serum (IgG). ABIN1 was subsequently detected by Western Blot using the anti-ABIN1 antibody (Cat# 68-0001-100). In order to show that all ABIN1 was immunoprecipitated from the input cell extract, a Western Blot was carried out using anti-ABIN1 antibody (Cat# 68-0001-100) on 20 µg of the cell extract following immunoprecipitation - 'Supernatant' - and no ABIN1 could be detected. This demonstrates that 0.3 µg of the anti-ABIN1 antibody (Cat# 68-0001-100) can completely deplete ABIN1 from 1 mg of Raw 264.7 total cell extract.



www.ubiquigent.com
Dundee, Scotland, UK

ORDERS / SALES SUPPORT

International: +1-617-245-0020
US Toll-Free: 1-888-4E1E2E3 (1-888-431-3233)
Email: sales.support@ubiquigent.com

UK HQ and TECHNICAL SUPPORT


International: +44 (0) 1382 381147 (9AM-5PM UTC)
US/Canada: +1-617-245-0020 (9AM-5PM UTC)
Email: tech.support@ubiquigent.com

Email services@ubiquigent.com for enquiries regarding compound profiling and/or custom assay development services.

© Ubiquigent 2014. Unless otherwise noted, Ubiquigent, Ubiquigent logo and all other trademarks are the property of Ubiquigent, Ltd.

Limited Terms of Use: For research use only. Not for use in humans or for diagnostics. Not for distribution or resale in any form, modification or derivative OR for use in providing services to a third party (e.g. screening or profiling) without the written permission of Ubiquigent, Ltd.

Lot-specific COA version tracker: v1.0.0

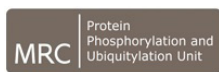

ABIN1 (mouse; full length), pAb

Alternate Names: A20-binding inhibitor of NF-kappa-B activation 1, HIV-1 Nef-interacting protein, Nef-associated factor 1, Nip40-1, Virion-associated nuclear shuttling protein, TNFAIP3-interacting protein 1

Cat. No. 68-0001-100
Lot. No. 30238

Quantity: 100 µg
Storage: -20°C

FOR RESEARCH USE ONLY

NOT FOR USE IN HUMANS

CERTIFICATE OF ANALYSIS

Page 2 of 2

Background

Continued from page 1

General References:

Gao L, Coope H, Grant S, Ma A, Ley S and Harhaj E (2011) ABIN1 Protein Cooperates with TAX1BP1 and A20 Proteins to Inhibit Antiviral Signaling. *J Biol Chem* **286**, 36592–36602.

Nanda S, Venigalla R, Ordureau A, Patterson-Kane J, Powell D, Toth R, Arthur J, and Cohen P (2011) Polyubiquitin binding to ABIN1 is required to prevent autoimmunity. *J Exp Med* **208**, 1215–1228.

Verstrepena L, Carpentiera I, Verhelsta K and Beyaert R (2009) ABINs: A20 binding inhibitors of NF-κB and apoptosis signaling. *Biochem Pharmacol* **78**, 105–114.

Zhou J, Wu R, High A, Slaughter C, Finkelstein D, Rehg J, Redecke V and Häckera H (2011) A20-binding inhibitor of NF-κB (ABIN1) controls Toll-like receptor-mediated CCAAT/enhancer-binding protein β activation and protects from inflammatory disease. *Proc Natl Acad Sci USA* **108**, 998–1006.

Application Reference:

Nanda S, Venigalla R, Ordureau A, Patterson-Kane J, Powell D, Toth R, Arthur S, Cohen P (2011) Polyubiquitin binding to ABIN1 is required to prevent autoimmunity. *J Exp Med* **208**, 1215-28.



www.ubiquigent.com
Dundee, Scotland, UK

ORDERS / SALES SUPPORT

International: +1-617-245-0020
US Toll-Free: 1-888-4E1E2E3 (1-888-431-3233)
Email: sales.support@ubiquigent.com

UK HQ and TECHNICAL SUPPORT

International: +44 (0) 1382 381147 (9AM-5PM UTC)
US/Canada: +1-617-245-0020 (9AM-5PM UTC)
Email: tech.support@ubiquigent.com

Email services@ubiquigent.com for enquiries regarding compound profiling and/or custom assay development services.

© Ubiquigent 2014. Unless otherwise noted, Ubiquigent, Ubiquigent logo and all other trademarks are the property of Ubiquigent, Ltd.

Limited Terms of Use: For research use only. Not for use in humans or for diagnostics. Not for distribution or resale in any form, modification or derivative OR for use in providing services to a third party (e.g. screening or profiling) without the written permission of Ubiquigent, Ltd.

Lot-specific COA version tracker: v1.0.0